

DrägerService®

Dräger

Field Service Procedure

Part Number: SP00258

Rev: -

Date: 18 Dec 2002

© 2002 Draeger Medical, Inc.

Installation of APL Valve on Fabius GS

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

Installation of APL Valve on Fabius GS

Replacement Procedure

Preparing the Breathing System:

1. Set machine power switch to the OFF (⦿) position and close all flow control valves.
2. Remove APL valve from the breathing system.
3. Record the serial number of the original APL Valve and replacement APL Valve for subsequent documentation step. See Figure 1 for serial number locations.
4. Install new APL Valve (P/N MK00625) securely into place with the retaining nut so that the pressure setting indicator is facing the Fabius GS operator. See Figure 1.
5. Locate product Operator's Manual and insert APL Valve Supplement behind current documents in binder assembly.

6. Discard original APL Valve.

7. If performing a software update at time of APL Valve replacement, refer to Software Upgrade Service Procedure for performance tests required. Otherwise, proceed to step 8.

Performance Tests:

8. Perform the following steps outlined in **Service Procedure SP00225 (REV. F)**.

Sections:

7.9 Breathing System

7.21 Ventilator

9. Complete a corresponding Service Report indicating serial number of APL Valves removed and installed, and test results (pass/fail).

Replacement Procedure (continued)

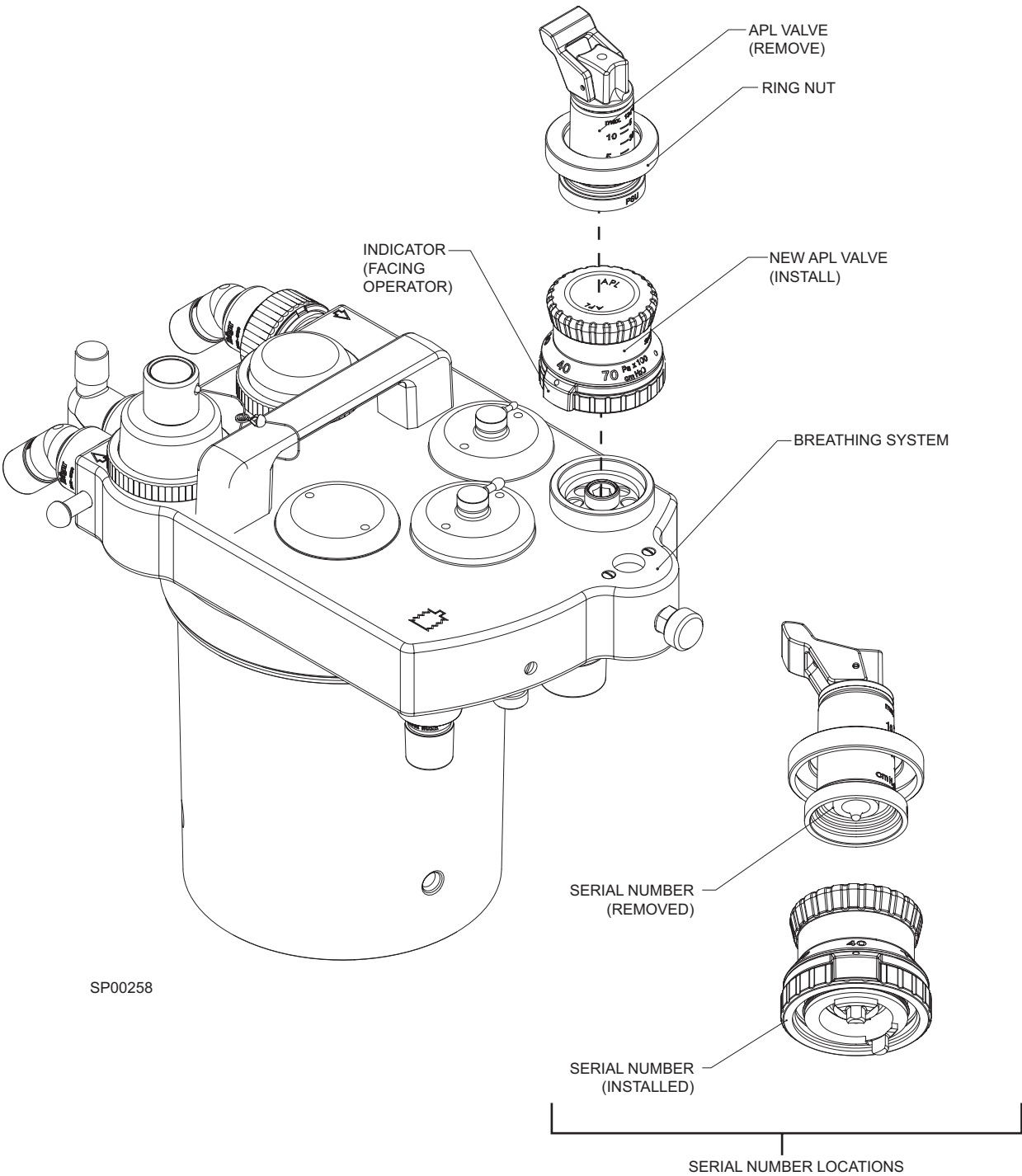


FIGURE 1. APL Valve Replacement

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

DrägerService®

Dräger

DrägerService is a division of
Draeger Medical, Inc.
3122 Commerce Drive
Telford, PA 18969
Tel: (215) 721-5402
(800) 543-5047
Fax: (215) 721-5784
Web: www.draegermedical.com
Printed in the U.S.A.